

Best Stairlift Brands in Singapore: Models, Fit and Support

By [DirectHome Team](#)

Updated 16 September 2026

8 min read

[Explore our stairlifts service ↗](#)


For a straight indoor stair, compare models such as Acorn 130, Stannah Siena and Bruno Elan. For a turning route, Acorn 180, Stannah Siena curved, Handicare Freecurve and Otolift Modul-Air Smart offer different rail and assistance arrangements. The best choice depends on the rider's weight, ability to sit and transfer, and the measured staircase; the Singapore installation and service arrangement must support that exact configuration.

Need advice for your home?

[WhatsApp our team ↗](#)



A model and its local support go together

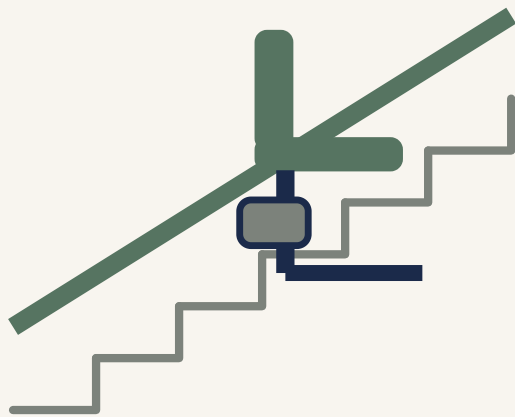


Exact model

Local support
Written scope

Schematic • not to scale

A stairlift proposal needs both an exact model and a local support arrangement. A familiar brand name does not establish either scope.

Why Brand Actually Matters for a Stairlift

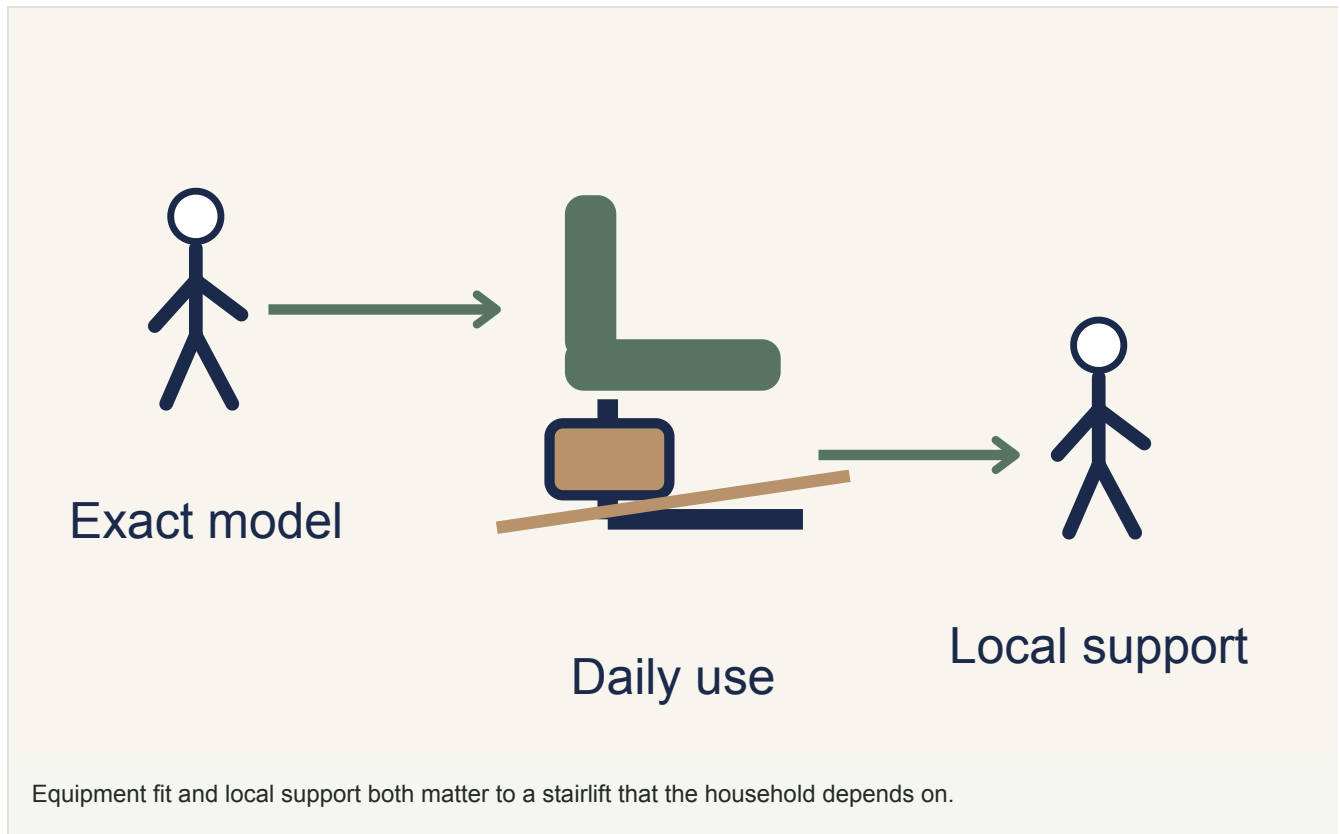

The chair and the local service arrangement are two separate purchases. A suitable seat, controls and rail determine whether the user can complete the journey. The supplier's support determines what happens when the equipment stops working or a part is needed. Compare both before selecting a brand.

Start with a short list of models that fit the actual staircase and user. Ask each supplier to name the Singapore installation and maintenance contractors, show the complete configuration and provide written support terms. A familiar badge cannot substitute for those answers.

Manufacturer reference configurations. Confirm the rating, options and support of the actual Singapore supply; these are not measured comfort rankings or assurances of independent transfer.


Model to compare	Documented distinction	Household decision it helps resolve
Acorn 130	Straight; 127 kg standard or 159 kg HD; powered-swivel configuration offered	A heavier rider needs the correct version, not a brand-wide capacity
Stannah Siena	Published reference capacities: 158 kg straight, 136 kg curved; assistance depends on configuration	The same chair name has different ratings on straight and turning routes
Bruno Elan	Straight indoor; 300 lb (about 136 kg); optional power swivel and footrest	A straight route with specific folding or turning assistance needs

Model to compare	Documented distinction	Household decision it helps resolve
Acorn 180	Curved; 120 kg; modular single rail	A turning route, provided the rider and local configuration meet the rating
Handicare Freecurve	Custom single rail; optional Active Seat stand assistance and powered swivel	A rider who needs help with standing as well as travelling upstairs
Otolift Modul-Air Smart	Modular rail; automatic seat swivel and automatic footrest folding described	A turning route where manual folding or turning is difficult


For a rider weighing 130 kg, the Acorn 180's 120 kg rating rules out that listed configuration. The Siena curved reference of 136 kg is a different starting point, subject to the supplied model, seat fit and an assessment of the rider's transfer. Conversely, a suitable straight staircase opens up the Acorn 130 HD and other straight models. Capacity narrows the choice before upholstery or promotional prices matter.

Consider the space left for other stair users as well as whether the chair can fit. [stairlift options for narrow stairs →](#)

Stannah: Compare the Seat and Operating Options

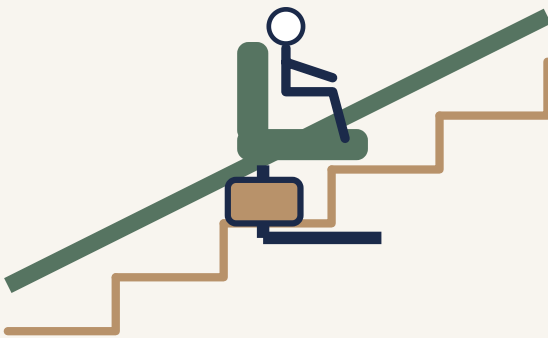

Start + stop

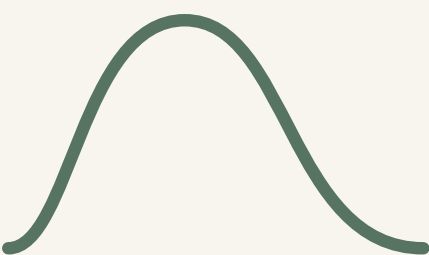

Ride character

Try the selected model

An actual demonstration helps the user judge the start, travel and stop; this is not a measured performance graph.

For a Stannah proposal, identify the chair and rail model rather than comparing the brand in the abstract. The [260 Siena manual](#), for example, describes specific seat, swivel and operating arrangements. Other configurations can differ, so bring the intended user's needs to the demonstration.

The published Siena reference is rated to 158 kg on a straight stair and 136 kg on a curved route; do not carry the higher figure across both. [Powered Easy-Turn](#) and [powered footrest folding](#) are configuration options. They can reduce the turning or reaching required, but the user still needs a suitable seat and transfer arrangement at each landing.


Try the belt, controls, footrest and landing transfer with suitable assistance. Determine which powered options are included and whether they solve a real difficulty. Price those options in the local quotation; an overseas specification or warranty does not establish the terms available for your home.

Want an exact price for your home?

Ranges only go so far. WhatsApp us a photo of your space and we'll quote your actual job. No obligation.

[WhatsApp us](#) [Try the free cost estimator](#)

Acorn: Distinguish the Straight and Curved Model

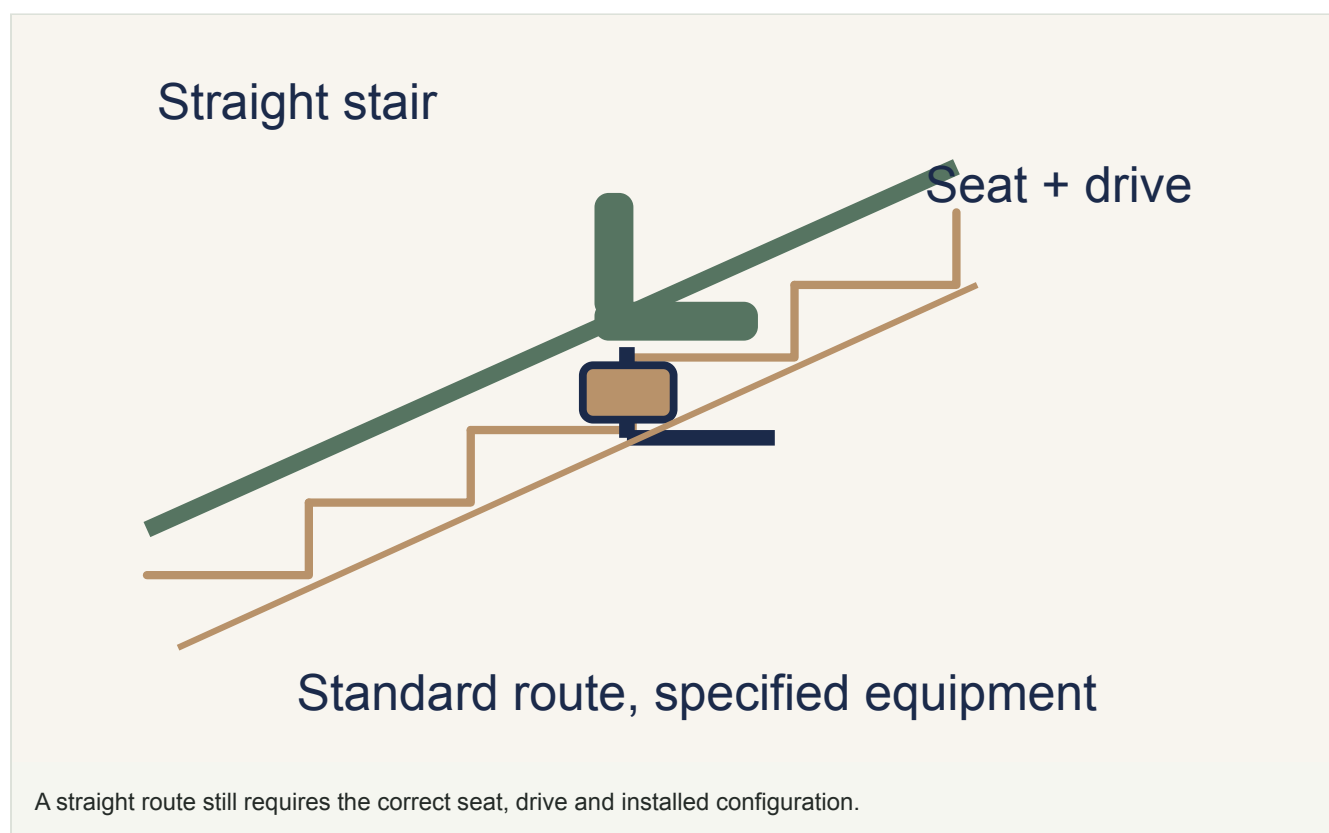


For an Acorn offer, check whether it is for a straight or curved route and record the exact model and options. The manufacturer's model-specific capacity information illustrates why a single brand-wide rating would be misleading. Use the documentation for the locally supplied configuration.

The Singapore Acorn pages make the distinction concrete: the 130 straight model lists 127 kg standard and 159 kg in the HD version, while the 180 curved model lists 120 kg. The 180 uses a modular single rail. Modularity describes how the route is assembled; it does not by itself promise faster Singapore delivery, a cheaper installation or a fit on every bend.

Compare the complete rail route, landing arrangement and user trial with the other offers. Ask how the chair is charged and how a fault is reported. Price and speed of supply are useful only after suitability, installation and ongoing support have been established.

Separate the equipment, installation and ongoing costs when comparing offers. [stairlift costs in Singapore](#) →

Bruno: Match the Product to Its Environment

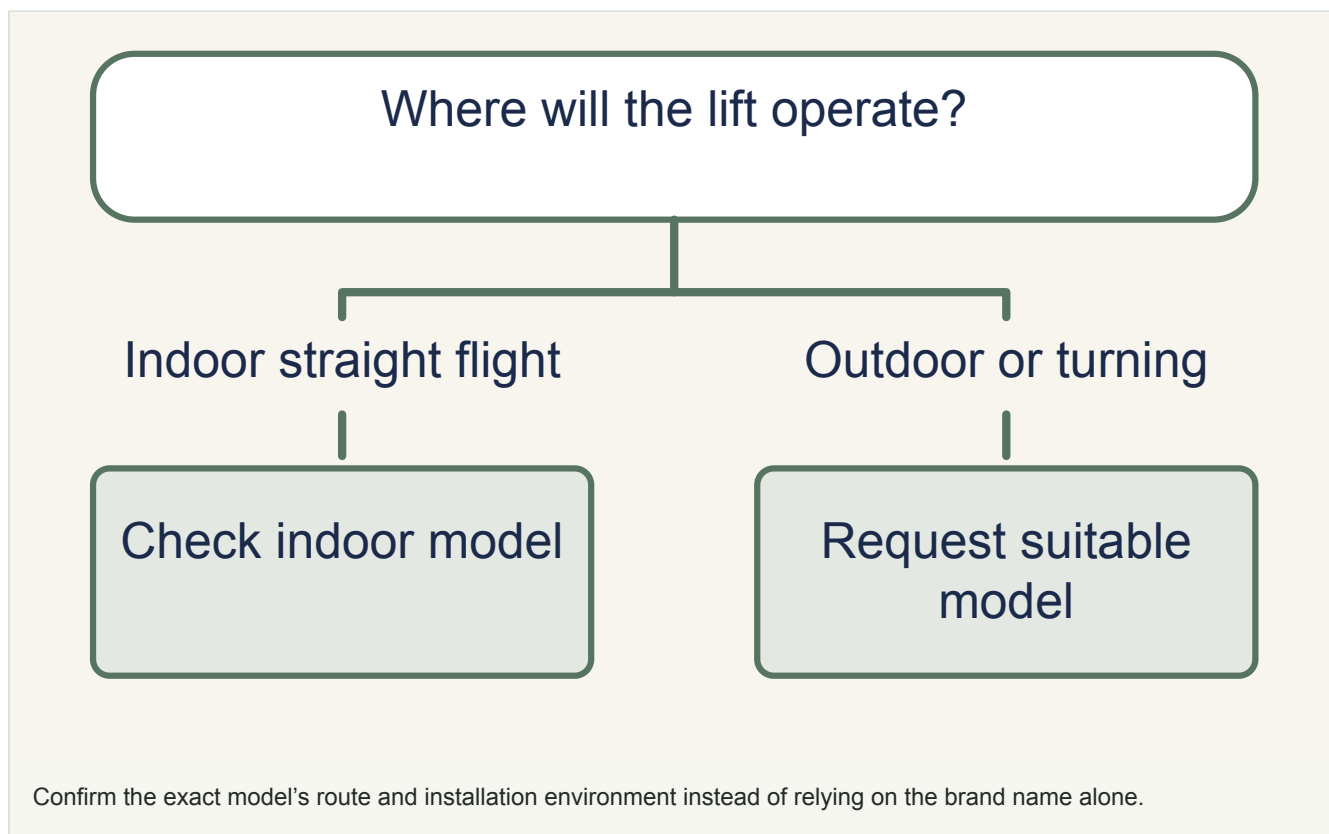

The Bruno Elan is a straight indoor model. A staircase with turns or an exposed outdoor location needs a different suitable product. Confirm the exact rail, seat and environment before comparing the quotation with another manufacturer's chair.

The Elan's 300 lb rating is approximately 136 kg. Its optional power swivel and power footrest are useful to trial if turning the seat or reaching down is difficult. This is a specific straight indoor alternative; an outdoor staircase needs an appropriately rated outdoor product, not an indoor Elan

protected only by a cover.

Ask the local provider to identify who holds the installation and maintenance responsibility, which parts are supported and the written response arrangements. A manufacturer's overseas presence does not establish local stock or attendance times. Keep those commercial commitments alongside the surveyed layout.

Handicare and Otolift: Compare the Complete Curved Route


A brand shortlist starts with the stair



Straight



Turning

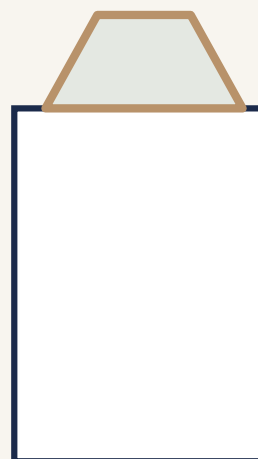

Outdoor

Match the exact model to the route

Schematic • not to scale

A straight route, a turning route and an outdoor setting create different equipment questions. Confirm the appropriate model before comparing brands.


If a supplier proposes a Handicare or Otolift system for a curved staircase, compare its occupied swept path with the alternatives. A narrow rail profile does not by itself establish knee clearance, headroom or space for a caregiver at the landing. Have those dimensions shown on the same surveyed staircase.

The [Handicare Freecurve](#) uses a custom single rail and offers an Active Seat option to assist standing up and sitting down, as well as powered swivel and footrest options. The [Otolift Modul-Air Smart](#) uses a modular rail and describes automatic seat swivel and footrest folding. These are different kinds of help: a powered footrest removes a bending task, while stand assistance addresses part of getting out of the seat. Have the rider try the actual equipment with suitable assistance before treating either as the answer to a transfer difficulty.


For a tight turning stair, compare how the occupied chair travels and where it parks on the same measured route. Rail diameter alone does not tell you where knees, feet or a caregiver will be. If standing is the main difficulty, a Freecurve Active Seat trial may be more informative than a narrow-rail claim; if folding and turning are the difficult movements, trial the appropriate powered configurations across the shortlist.

Pay attention to where the chair parks and how the user gets off. Different rail arrangements can place the seat in quite different positions after a turn. Record the exact model, options, local support and installed scope rather than assuming all curved systems from a particular brand suit every tight bend.

Turn the Brand Shortlist Into Comparable Proposals

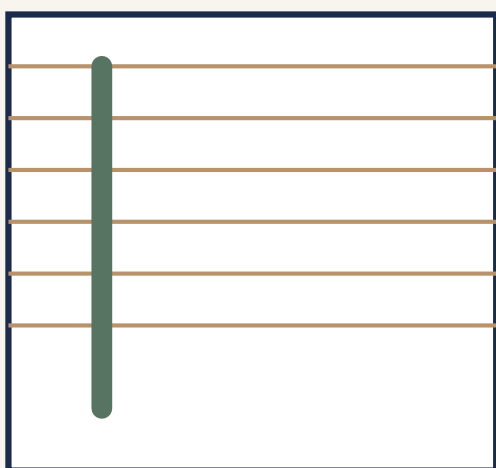

Straight flight

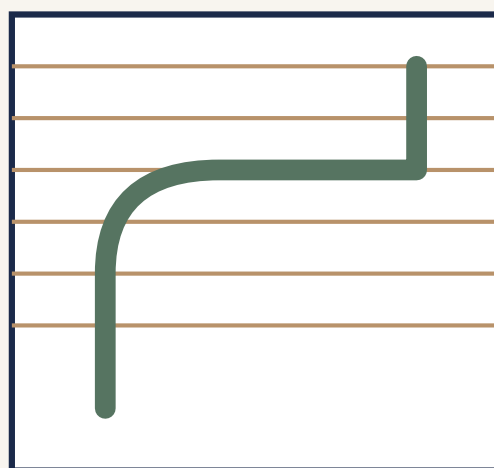

Turning staircase

Start brand shortlisting from the actual stair route before comparing suitable model offers.

Give each supplier the same user brief and staircase information, then request a site assessment. Compare the models that pass that first suitability check through seat fit, transfer sequence, rail and remaining walking space, charging, support and complete cost. Choose the configuration that resolves those requirements for this household.

In Singapore, confirm [BCA registration for installation and maintenance and a valid PTO before use](#). If enquiring through [DirectHome](#), describe the access problem and request the responsible provider's written scope. Do not assume that an enquiry includes project supervision, technical certification or a particular manufacturer's quotation.

Frequently asked questions

Which stairlift brand is the cheapest in Singapore?

For a suitable straight staircase, the [Acorn 130 is a starting point for a budget comparison](#): it is described as the least expensive model in that range. That does not prove Acorn is cheaper than every other brand. Compare complete local offers for the same route, required options and support; a straight model cannot replace a continuous curved route where the rider cannot transfer safely.

Which stairlift brand is best for a curved staircase?

Compare actual curved configurations: Acorn 180 has a 120 kg rating and modular rail; the listed Stannah Siena curved reference is 136 kg; Handicare Freecurve offers custom rail and optional stand assistance; Otolift Modul-Air Smart uses modular rail with automatic folding and swivel features. Surveyed clearance, the supplied model and the rider's transfer determine which can work.

Why does the brand matter so much for a stairlift?

The model affects fit and operation; the local provider affects installation, parts and service. Both need to be evaluated separately.

Is Stannah worth paying more for?

An extra cost can be worthwhile if the specified Siena seat and powered options solve a difficulty the rider has, or if its rated configuration fits where another does not. Trial those differences and compare written Singapore support terms. The Stannah badge alone does not establish better comfort or value.

More on stairlifts

[Article](#) [Can Two Family Members Use the Same Stairlift?](#)

[Article](#) [Perch Stairlifts: An Option When Sitting Is Difficult?](#)

[Article](#) [Stairlifts in HDB Maisonettes: Approval and Installation](#)

[Related service](#) [Home Lifts](#)

Our service

Stairlifts for your home

Explore the service, then discuss what your home needs with our team.

[View service details ↗](#)

Still weighing your options?

Plan with the buying guide

Bring the choices, costs and questions together before deciding what to do next.

[Read the stairlifts buying guide →](#)

[More stairlifts articles →](#)

Ready to discuss your home? [WhatsApp our team](#) or [request a quote](#).